

ARROW®


Critical Care Products

Arterial Catheterization Systems

Featuring
Spring-Wire
Guide Placement
Technique

Positive placements are easier with the first catheter designed specifically for arterial use.

- The Arrow Radial Artery Catheter features a patented integral spring-wire guide technique.
- An extra-short bevel on the introducer needle and the precision fit between the needle bevel heel and catheter tip facilitate accurate placement within the vessel lumen.

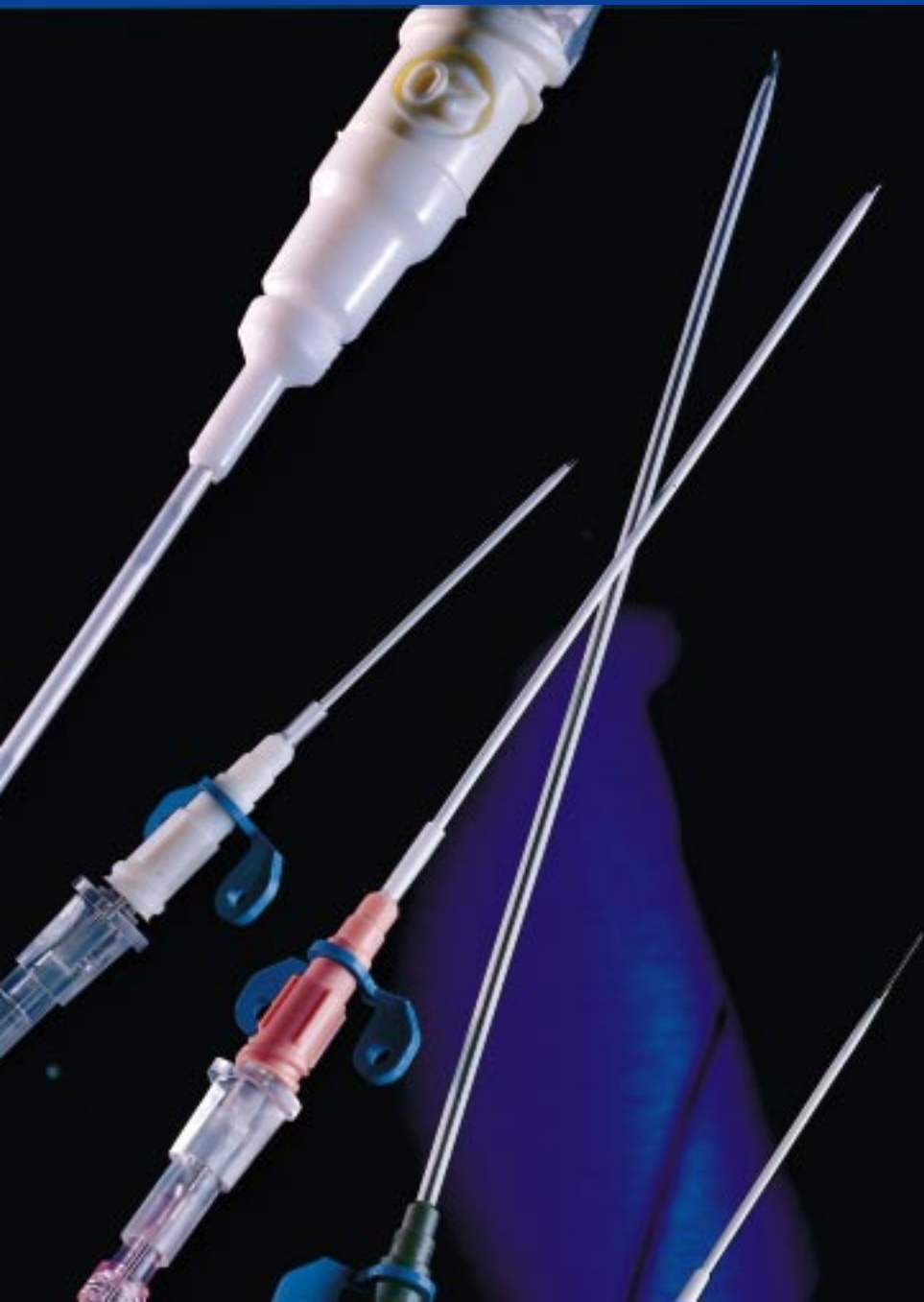

Arterial Catheterization Systems

Save time when performing radial, femoral and brachial artery catheter placements.

Integral Wire Products (Integral wire systems inhibit blood exposure)

Product No.	Description	CE	kits/ case
AK-04018††	Arrow 18 Ga. x 1-3/4" (4.40cm) Radial Artery Catheterization Kit		10
AK-04020†	Arrow 20 Ga. x 1-3/4" (4.40cm) Radial Artery Catheterization Kit		10
AK-04020-C†	Arrow 20 Ga. x 1-3/4" (4.40cm) Radial Artery Catheterization Kit		10
AK-04120*†	Arrow 20 Ga. x 1-1/2" (3.81cm) Radial Artery Catheterization Kit		10
AK-04122*†P	Arrow 22 Ga. x 1-3/8" (3.49cm) Radial Artery Catheterization Kit		10

Product No.	Description	CE	kits/ case
FA-04014†	Arrow 14 Ga. x 6" (16cm) Arterial/Venous Catheterization Set	CE	25
FA-04016†	Arrow 16 Ga. x 6" (16cm) Arterial/Venous Catheterization Set	CE	25
FA-04018††	Arrow 18 Ga. x 4-1/4" (10.8cm) Arterial Catheterization Set	CE	25
FA-04020††	Arrow 20 Ga. x 4-1/4" (10.8cm) Arterial Catheterization Set	CE	25
RA-04018††	Arrow 18 Ga. x 1-3/4" (4.4cm) Radial Artery Catheterization Set	CE	50
RA-04020†	Arrow 20 Ga. x 1-3/4" (4.4cm) Radial Artery Catheterization Set	CE	50
RA-04120*†	Arrow 20 Ga. x 1-1/2" (3.81cm) Radial Artery Catheterization Set	CE	50
RA-04122*†P	Arrow 22 Ga. x 1-3/8" (3.49cm) Radial Artery Catheterization Set	CE	50

QuickFlash® Products		CE	units/ case
AK-04220*†	Arrow 20 Ga. x 1-1/2" (3.81cm) QuickFlash® Radial Artery Cath. Kit		10
RA-04220†	Arrow 20 Ga. x 1-1/2" (3.81cm) QuickFlash® Radial Artery Cath. Set	CE	50
RA-04220-W†	Arrow 20 Ga. x 1-1/2" (3.81cm) QuickFlash® Radial Artery Cath. Set	CE	50

Safety Products		CE	units/ case
AK-04020-SPC	Arrow 20 Ga. x 1-3/4" (4.4cm) Radial Artery Catheterization Kit		10
RA-04020-SP	Arrow 20 Ga. x 1-3/4" (4.4cm) Radial Artery Catheterization Set		50

All-in-one design.


- Blood flashback in the clear hub of the introducer needle shows successful entry.
- Simple actuating lever advances the spring-wire guide through the needle lumen into the vessel.

Excellent stability and indwelling characteristics.

- Arrow's polyurethane construction softens *in situ*, offers excellent thromboresistance and

biocompatibility, and minimizes trauma to vessel walls.

- The catheter is reinforced at the hub juncture to minimize kinking at this maximum stress point.

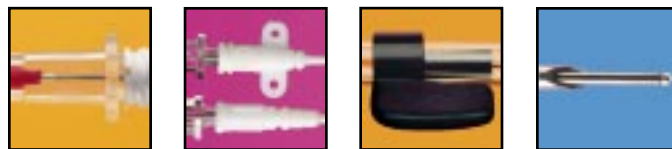

Traditional Seldinger Products

Product No.	Description	CE	kits/ case
AK-04510††	Arrow 20 Ga. x 5" (12.7cm) Arterial Line Kit		10
NS-04100††	Arrow 20 Ga. x 5" (12.7cm) Arterial Line Kit		10
UM-04018*†	Arrow 18 Ga. x 6" (16cm) Arterial Line Kit		10

Product No.	Description	CE	sets/ case
AV-04020††	Arrow 20 Ga. x 1-3/4" (4.4cm) Arterial/Venous Catheterization Set	CE	100
GH-04120*†	Arrow 20 Ga. x 2-1/2" (6.35cm) Radial Artery Catheterization Set	CE	25
GH-04120-E*†	Arrow 20 Ga. x 3-1/16" (7.78cm) Radial Artery Catheterization Set	CE	25
GH-04122*†	Arrow 22 Ga. x 1-3/8" (6.35cm) Radial Artery Catheterization Set	CE	25
GH-04125*	Arrow 20 Ga. x 3-1/16" (7.78cm) Radial Artery Catheterization Set	CE	25

Spring-Wire Guides

		CE	wires/ case
RW-04018	Arrow Positive Placement .018" (.45mm) x 5-1/4" (13.3cm) Spring-Wire Guide	CE	25
RW-04025	Arrow Positive Placement .025" (.64mm) x 5-1/4" (13.3cm) Spring-Wire Guide	CE	25

* Integral Suture Wings † Polyurethane †† FEP/Teflon® P Pediatric Size

Caution: U.S. Federal Law limits this device to sale by or on order of a physician. Contents of unopened, undamaged package are sterile. Disposable. Refer to package insert for current warnings, indications, contraindications, precautions, and instructions for use.

Arrow Radial Artery Catheter, US Patent No. 5,246,426.

Contact Arrow International, Inc. for reference information on the above.

To locate an Arrow Distributor or for ordering information contact customer service at:

United States

Tel: 800.523.8446

Fax: 800.343.2935

International

Tel: +1.610.320.3931

Fax: +1.610.320.3932

Distribution Worldwide:

Arrow Sales offices are located in Canada, Czech Republic, France, Germany, Greece, Holland, Japan, Mexico, Slovakia, South Africa, Spain, and the United States.

For complete contact information on Arrow International, Inc. offices, visit our website at www.arrowintl.com.

Corporate Headquarters address: Arrow International, Inc., 2400 Bernville Road, Reading, PA 19605 Tel: 610-378-0131

