


Bipolar Pacing Catheters and Pacing Kits

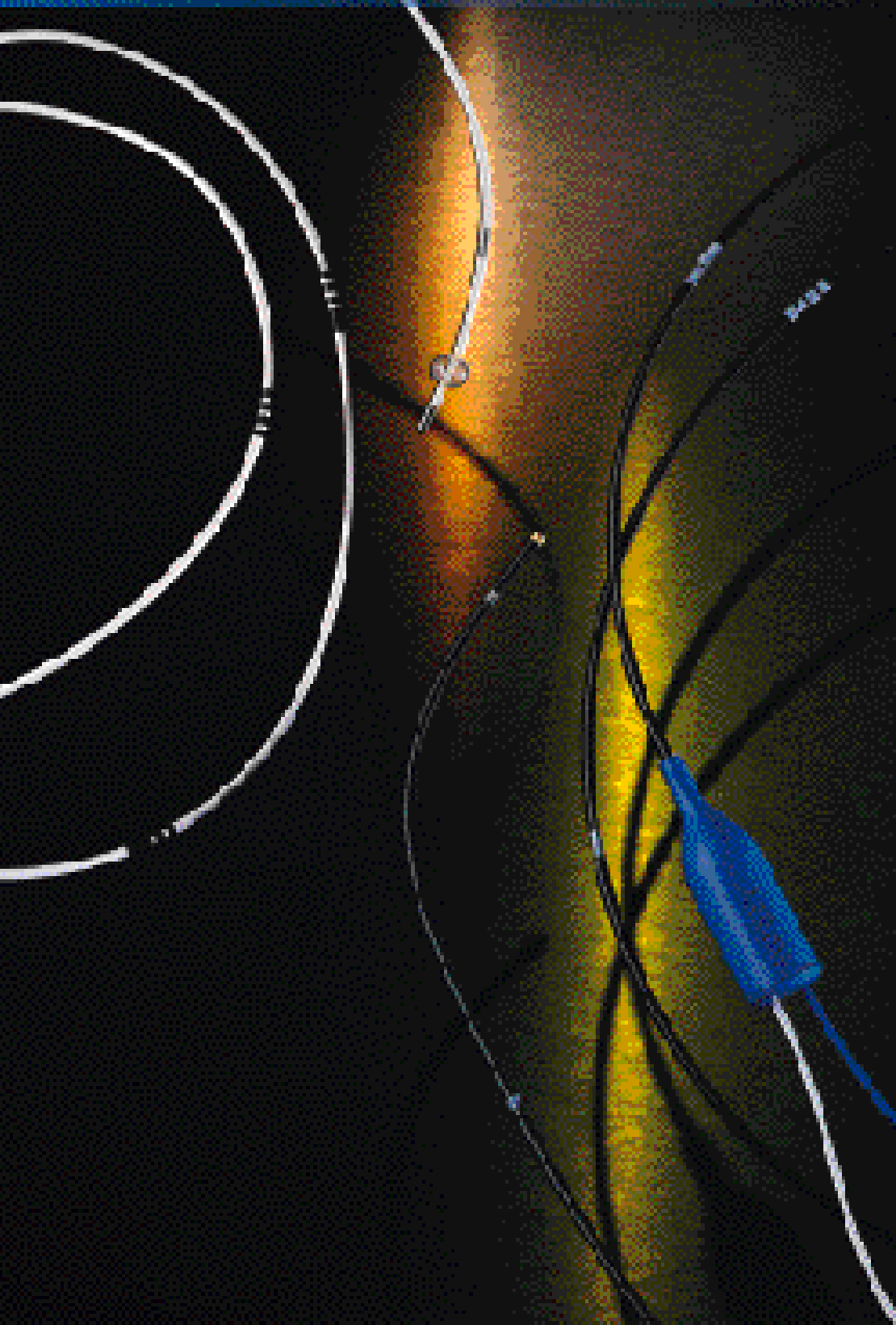


A
B
P
O
W
A

A complete line of
pacing catheters and
insertion kits for all
of your pacing
requirements.


Arrow Sheath Introducer Kits bring added convenience to pacing procedures.

When pacing the heart, having the right tools at hand is essential. Arrow, the world leader in sheath introducers, offers a complete line of bipolar pacing catheterization kits designed to help make the placement of pacing catheters easier and more convenient than ever.

Each kit includes one of Arrow's hardworking Bipolar Pacing Catheters: our Torque Directed Catheter, designed for superior maneuverability and mechanical stability, or our Balloon Flotation Catheter, which floats into the correct position at bedside without fluoroscopy.

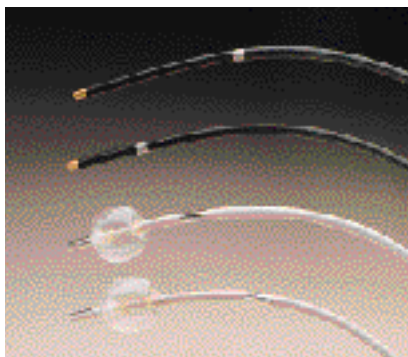

Bipolar Catheter Tips

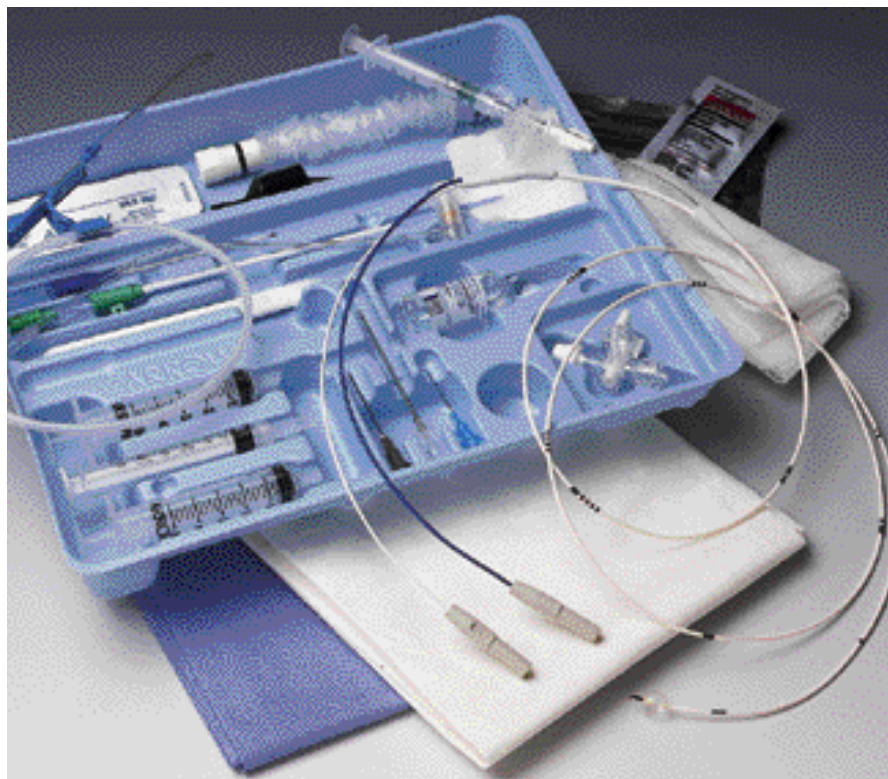

Balloon Flotation Catheter with Introducer Kit

All kit components are matched to fit.

Everything you need comes packaged in a single, compact kit. A variety of appropriate needles, syringes, prep solution, lidocaine and drapes are included, along with the catheter, a sheath dilator assembly, a side port/hemostasis valve and additional items. All PSI components have been matched to fit the Arrow Bipolar Pacing Catheter.

Easier to insert at the patient's bedside.

Because Arrow offers a Balloon Catheterization Kit, bipolar pacing procedures can be performed conveniently in the patient's room or in the ICU or CCU without fluoroscopy.¹ This is especially important when pacing needs to be established quickly and efficiently.

¹ Final position of pacing catheter should be confirmed by x-ray.



Arrow Sheath Introducer

Consider downsizing.

Arrow Torque Directed Catheters offer superior torqueability. Try our 5 Fr. compared to a competitor's 6 Fr. catheter. With the Arrow Torque Directed Catheter, downsizing does not compromise your torqueability.


Each kit contains the following components:

- One: Bipolar Pacing Catheter
- One: Sheath/Dilator Assembly
- One: Side Port/Hemostasis Valve with Locking Mechanism
- One: 3-Way Luer-Lock Stopcock
- One: 035" (45cm) "J" Tip Wire Guide
- One: Cath-Gard® with Sterile Tape
- One: 1.25cc Luer-Lock Syringe
- One: 3cc Luer-Lock Syringe
- One: 5cc Luer-Lock Syringe
- One: 25 Ga. 1" Needle
- One: 22 Ga. 1-1/2" Needle
- One: 18 Ga. 2-1/2" T.W. Needle
- One: 18 Ga. 2-1/2" Catheter over Needle
- One: Disposable Scalpel (#11 Blade)
- One: Straight Silk Suture
- Three: Povidone-Iodine Swabsticks*
- One: Povidone-Iodine Ointment* 1gm 10%
- One: Xylocaine 5 ml 1%*
- Two: 2" x 2" Gauze Pads
- Four: 4" x 4" Gauze Pads Assembly
- One: ECG Adapter
- One: Fenestrated Drape
- One: Non-Fenestrated Drape
- One: CSR Wrap

* For kits without medication items, order "-IK" product number. Available for export sale only.

Effective May 9, 2000 FDA regulation #21 CFR 898 restricts sales in the U.S. to Bipolar Catheters only with shrouded (covered) connector pins.

Complete Procedural Kits†

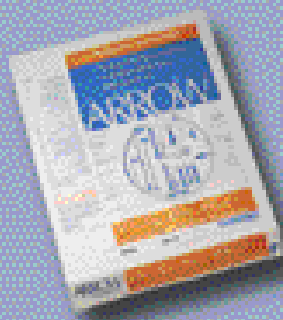

Includes Bipolar Pacing Catheter and Sheath Introducer

Product No.	French Size	Introducer Included	Recommended Useful Length	Inflated Balloon Capacity	Maximum Inflation Diameter	Electrode Width	Electrode Spacing	Medication in Kit
Balloon Flotation Catheter Kits								
AI-07155-KS	5	6 Fr.	110cm	3/4 cc	8mm	5mm	25mm	yes
AI-07155-IK*	5	6 Fr.	110cm	3/4 cc	8mm	5mm	25mm	no
Torque Directed Catheter Kits								
AI-06210-KS	6	7 Fr.	110cm	n/a	n/a	2mm	10mm	yes
AI-06210-IK*	6	7 Fr.	110cm	n/a	n/a	2mm	10mm	no

Bipolar Pacing Catheters†



Catheter available as separate item

Product No.	French Size	Recommended Introducer Size	Recommended Useful Length	Inflated Balloon Capacity	Maximum Inflation Diameter	Electrode Width	Electrode Spacing
Balloon Flotation Catheters							
AI-07154*	5	6 Fr.	105cm	3/4 cc	8mm	5mm	15mm
AI-07154-S	5	6 Fr.	105cm	3/4 cc	8mm	5mm	15mm
AI-07155*	5	6 Fr.	110cm	3/4 cc	8mm	5mm	25mm
AI-07155-S	5	6 Fr.	110cm	3/4 cc	8mm	5mm	25mm
HANDS-OFF® Balloon Flotation Catheter							
AH-07155*	5	6 Fr.	110cm	3/4 cc	8mm	5mm	25mm
Torque Directed Catheters							
AI-05210*	5	5 Fr.	110cm	n/a	n/a	2mm	10mm
AI-05210-S	5	5 Fr.	110cm	n/a	n/a	2mm	10mm
AI-05225*	5	5 Fr.	110cm	n/a	n/a	2mm	25mm
AI-06210*	6	6 Fr.	110cm	n/a	n/a	2mm	10mm
AI-06210-S	6	6 Fr.	110cm	n/a	n/a	2mm	10mm
AI-06225*	6	6 Fr.	110cm	n/a	n/a	2mm	25mm
Balloon Flotation Pacing Catheter with Two Pressure Lumens¹							
AI-17210*	7	7 Fr.	110cm	1-1/2cc	12mm	5mm	10mm
AI-17210-S	7	7 Fr.	110cm	1-1/2cc	12mm	5mm	10mm

NOTE: Products with "S" or "KS" have shrouded connector pins for domestic sale.

† All catheter tips have right heart curves.

* Indicates product for export sale only.

¹ Electrodes at 14 and 15cm from catheter tip with ports at distal and 29cm.

Arrow Bipolar
Pacing Catheters
and Pacing Kits

Arrow also offers separately packaged individual Percutaneous Sheath Introducer Kits and Sets for use with Arrow Bipolar Pacing Catheters and other 5, 6 and 7 Fr. catheters.

Product No.	Sheath French Size	Sheath Material	Spring Wire Guide	Type and Quantity/case
To introduce 5 Fr. Catheters				
AK-09601	6 Fr.	FEP	.035"	10 Kits/case
SI-09600	6 Fr.	FEP	.035"	10 Sets/case
To introduce 6 Fr. Catheters				
AK-09701	7 Fr.	FEP	.035"	10 Kits/case
SI-09700	7 Fr.	FEP	.035"	10 Sets/case
To introduce 7 Fr. Catheters				
AK-09800	8.5 Fr.	Poly	.035"	10 Kits/case
AK-09801	8.5 Fr.	Poly	.035"	10 Kits/case
SI-09808	8 Fr.	FEP	.035"	10 Sets/case

Arrow International has a complete line of single-lumen and multi-lumen products to meet your central venous access needs. Please contact your local Arrow Representative or contact us by calling 1-800-523-8446 or 610-378-0131. We invite you to visit us on the world wide web at www.arrowintl.com for valuable information about all of our product lines.

For information regarding distribution in your area, contact:



www.arrowintl.com

UNITED STATES

Arrow International, Inc.
2400 Bernville Road
Reading, PA 19605, USA
Phone: 610-378-0131
Toll-free: (800) 523-8446
Fax: 610-478-3199
Orders-only toll-free fax: (800) 343-2935
Email: Customer.Service@arrowintl.com

AFRICA

Arrow Africa (Pty) Ltd.
Cambridge Commercial Park
22 Witkoppen Road, Paulshof Extension
Sandton 2054, Republic of South Africa
Phone: (11) 807 4887
Fax: (11) 807 4994
Email: tlc@arrowafrica.co.za

CANADA

Arrow Medical Products, Ltd.
2300 Bristol Circle, Unit 1
Oakville, Ontario L6H 5S3
Phone: (905) 829-9473
Toll-free: (800) 387-7819
Fax: (905) 829-9414
Email: arrow.canada@arrowintl.com

EUROPE

Arrow Deutschland GmbH
Justus-von-Liebig-Strasse 2
D-85435 Erding, Germany
Phone: 08122 9820-0
Fax: 08122 40384
Email: service@arrow-deutschland.de

Arrow France S.A.
Atlantic Parc, « Les Pyramides » No. 11
Route de Pitoys, P.A. de Maignon
64600 Anglet, France
Phone: 05 59 31 34 90
Fax: 05 59 31 34 91
Email: arrow.sa@wanadoo.fr

Arrow Hellas A.E.E.
230 Kifissias Avenue
Halandri, 152 31
Athens, Greece
Phone: 6777717 / 6777311
Fax: 6777911
Email: arrowhel@hol.gr

Arrow Holland Medical Products B.V.
Flevolaan 9A
NL 1382 JX
Weesp, The Netherlands
Phone: 294 280 620
Fax: 294 412 300
Email: c_service@arrow-holland.nl

Arrow Iberia S.A.
Avenida de Valgrande 14
28108 Alcobendas
Madrid, Spain
Phone: 916 621 267
Fax: 916 619 756
Email: arrow@arrow-iberia.es

Arrow International CR, a.s.
Pražská 209
500 04 Hradec Králové
Czech Republic
Phone: 049 575 9111
Fax: 049 575 9222
Email: jitka_vodstrcilova@arrow.cz

Arrow Slovensko Piešťany s.r.o.
Valová 49
921 01 Piešťany, Slovakia
Phone/Fax: 0838/77 25428
Email: arrow@arrow.sk

JAPAN

Arrow Japan, Ltd.
Nagaoka Bldg.
3-16, Kita-Otsuka 3-chome
Toshima-ku, Tokyo 170-0004, Japan
Phone: 03-5974-1701
Fax: 03-5974-1845
Email: arrowjapan@msn.com

Arrow Japan, Ltd.
Osaka Branch
9F Shinosaka Uenotoyo Bldg.
7-4-17 Nishinakanoshima
Yodogawa-ku, Osaka 532-0011, Japan
Phone: 06-6390-4554
Fax: 06-6390-4544
Email: arrowjapan@msn.com

MEXICO

Arrow Internacional de Mexico S.A. de C.V.
Molier No. 128, Col. Polanco
11560 Mexico City, Mexico
Phone: 5281-2291
Fax: 5282-1359
Email: arrowmex@mail.internet.com.mx

Caution: U.S. Federal Law limits this device to sale by or on order of a physician. Contents of unopened, undamaged package are sterile. For single use only. Do not sterilize or reuse.

Refer to package insert for current warnings, precautions, and instructions for use.

Contact Arrow International, Inc. for reference information on the enclosed.

©2000 Arrow International, Inc. All rights reserved. Printed in U.S.A.

BPP-C 11/00 8M

