


# ARROW®

CRITICAL CARE PRODUCTS DIVISION

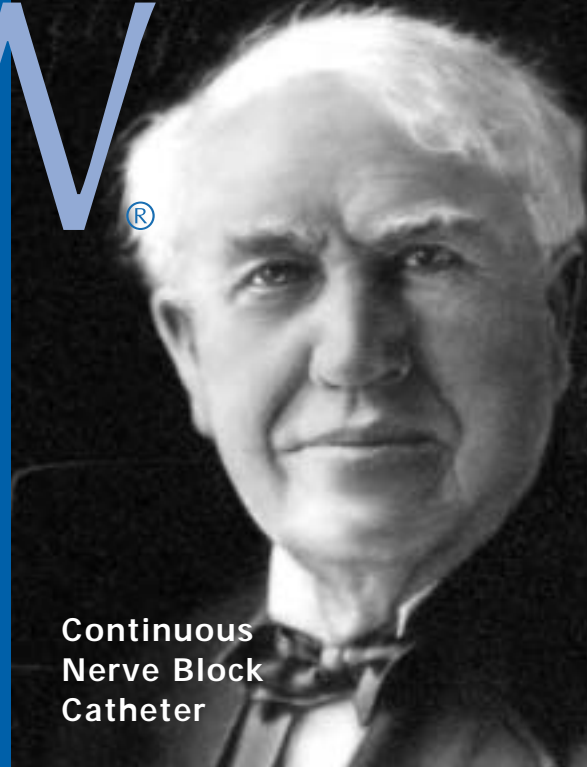

Continuous  
Nerve Block  
Catheter

## StimuCath™



**powered by innovation™**

Arrow salutes Thomas Edison—a shining example of how the power of innovation can improve the quality of lives around the globe.

# The Arrow StimuCath™

## ***A regional anesthesia product innovation***

- *Precise catheter placement*
- *Eliminates the likelihood of secondary block failure*
- *Easy to learn*
- *Simple to use*



# Continuous Nerve Block Catheter

*Arrow's unique StimuCath™ system allows nerve stimulation through both the needle and the catheter to assist in catheter placement for continuous plexus and peripheral nerve blocks.*



*Wire-reinforced technology.*



*Insulated needles ease insertion*



*StatLock®*



*SnapLock™*

## **Wire-reinforced technology enhances performance.**

- 19 Ga. bullet-tip catheter features Arrow's wire-reinforced catheter technology.
- Provides outstanding kink resistance and the ability to stimulate through the catheter.

## **Insulated needle eases insertion.**

- 17 Ga. blunt Tuohy tip insulated needle comes with catheter set.
- Provides improved tactile feel and catheter steerability.

## **StatLock® anchoring device secures adapter.**

- The StatLock® adhesive anchoring pad is custom designed for the Arrow SnapLock™ adapter.
- Securely affixes the catheter adapter comfortably to the patient.
- Decreases catheter displacements.

## **Catheter syringe adapter facilitates conduction.**

- StimuCath™ version of Arrow's popular SnapLock™ catheter/syringe adapter allows electric conduction to the catheter tip.
- Easy to use and decreases the risk of catheter disconnects.



*AB-05000 Set*

## StimuCath™ (Continuous Nerve Block Catheter Set)

Product No.	Description	Per Case	Product No.	Description	Per Case
AB-05000	<b>StimuCath™ Continuous Nerve Block Set</b> One: 19 Ga. StimuCath™ continuous nerve block catheter (90 cm) One: SnapLock™ adapter One: 17 Ga. x 3-1/2 in. (8.89 cm) Boezaart insulated plexus block TW needle with cm markings One: Connecting Cable One: Injection Needle: 18 Ga. x 1-1/2 in. (3.81 cm) One: 0.2 Micron Flat Filter One: Alcohol prep pad One: Dressing: StatLock** Anchoring Device 6.7 cm x 9.5 cm One: Tubing: 6 in. extension line with slide clamp	10	AB-05060-PK	<b>19 Ga. x 60 cm StimuCath™ Continuous Nerve Block Procedural Kit with Sharps Safety Feature</b> One: 19 Ga. StimuCath™ continuous nerve block catheter (60 cm) One: SnapLock™ adapter One: Injection Needle: 22 Ga. x 1-1/2" (3.81 cm) One: Injection Needle: 18 Ga. x 1-1/2" (3.81 cm) One: 17 Ga. x 3-1/2 in. (8.89 cm) insulated plexus block TW needle with cm marking One: Syringe: 10 cc Luer-Slip One: Syringe: 7 cc Luer-Slip Loss of Resistance One: Injection Needle: 25 Ga. x 5/8" (1.59 cm) and 3 cc Luer-Lock Syringe One: Alcohol prep pad One: Swabstick Foil Pack 23% Benzoin Tincture One: 5 ml Ampule HCl, 1% Lidocaine Solution One: 5 ml Ampule HCl, 1.5% Lidocaine w/Epinephrine 1:200,000 Solution One: Drape: 30" x 24" with 4" fenestration One: Towel: 13" x 18" One: CSR Wrap One: Dressing: StatLock** Anchoring Device 6.7 cm x 9.5 cm Two: Dressing: Tegaderm*** 10 cm x 12 cm One: Tubing: 6 in. extension line with slide clamp One: Connecting Cable One: Filter: 5 Micron Straw One: Filter: 0.2 Micron Flat Eight: Gauze Pad(s): 4" x 4" Three: Prep Sponge Swab(s) One: Prep Tray One: SharpsAway II*** Locking Disposal Cup	5
AB-05060	<b>StimuCath™ Continuous Nerve Block Set</b> One: 19 Ga. StimuCath™ continuous nerve block catheter (60 cm) Additional kit contents same as AB-05000 listed above.	10			
AB-00050	<b>Continuous Nerve Block Needle</b> One: 17 Ga. x 2 in. (5.00 cm) Boezaart insulated plexus TW needle with cm markings	10			
AB-00090	<b>Continuous Nerve Block Needle</b> One: 17 Ga. x 3-1/2 in. (8.89 cm) Boezaart insulated plexus TW needle with cm markings	10			
AB-00150	<b>Continuous Nerve Block Needle</b> One: 17 Ga. x 5-15/16 in. (15.0 cm) Boezaart insulated plexus TW needle with cm markings	10			

**Latex Free**

\* StatLock™ is a registered trademark of Venetec International Inc. Licensed under US Patent No. 5,192,274.

\*\* SharpsAwayII™: US Patent Nos. 6,123,193 and 6,276,527 B1.

\*\*\* Tegaderm™ is a registered trademark of 3M company.

Other Patents Pending.

Contact Arrow International, Inc. for reference information on the above.

**Caution:** U.S. Federal Law limits this device to sale by or on order of a physician. Contents of unopened, undamaged package are sterile. Disposable. Refer to package insert for current warnings, indications, contraindications, precautions, and instructions for use.

**Ordering Information (USA)**  
610.378.0131 or 800.523.8446  
FAX 800.343.2935 or 508.660.1819

**Ordering Information (International)**  
+1.610.320.3931  
FAX +1.610.320.3932

For information regarding distribution in your area, contact: Arrow International, Inc., 2400 Bernville Road, Reading, PA 19605, U.S.A.

**Distribution Worldwide:**

Arrow Sales offices are located in Czech Republic, France, Germany, Greece, Holland, Japan, Mexico, Portugal, Slovakia, South Africa, and Spain.

For complete contact information on Arrow International, Inc. offices, visit our website at [www.arrowintl.com](http://www.arrowintl.com).

To locate an Arrow distributor outside the U.S.A., contact International Customer Service at Tel: +1.610.320.3931 or FAX +1.610.320.3932.

